

STR

STR 2000 - STR 2000F



eryiğit

TIBBİ CİHAZLAR AŞ



ERYİĞİT Medical Device Company have ISO 9001, ISO 13485 and ISO 14001 quality management systems.



Our company, from product design, after-sales technical service at every stage of production till the "Quality Management" principles and guidelines as closely connected;

- Turnkey central sterilization snits (Autoclaves),
 - Steam sterilizers,
 - Operating aables and accessories and traction kits,
 - Operating ceiling lamp and patient consultation lamps,
 - Gynecological, urological and women on boards,
 - Hospital sterilization stainless steel equipment,
 - Life patient beds and transfer stretchers,
 - Otorinolarengoloji chair (Proctology tables),
- our products.



- EN 17665-1 Validation Test
- EN 285
- TIP Test certification
- ISO 9001
- ISO 13485
- 1984
- 1783



ERYİĞİT TABLES: General Surgery, Urology, Gynecology, Orthopedics, Ophthalmology, Neurosurgery, Ear, nose and throat, cardiology is designed for use in operations.

Operating tables and accessories

ERYİĞİT TABLES: General Surgery, Urology, Gynecology, Orthopedics, Ophthalmology, Neurosurgery, Ear, nose and throat, cardiology is designed for use in operations.

Table top tray, the head, the two-piece back, hip and leg portion including a two-part consists of six chapters. The table, up and down, the back part of Trendelenburg, reverse Trendelenburg, right - left hand control with the help of the lateral movements of the moving wire electro - hydraulic system.

Manually moving the head and foot part of the request may be extended up to 5 cm table height. Can be remove if you want.

Operating table (General purpose)

Foot section consists of two parts, 90 degree down, 180 degree left – right, 30 degree up and is able to move. The lower part of the plate and the moving body, covered with stainless steel telescopic panels are protected against external influences, and all types of liquids. These panels are removable . On the top plate, 5 cm thick, antistatic and antibacterial or 8 cm thick, viscoelastic polyurethane cushions are available. cushions are removable and easy clean. Suitable for use with C-arm X-ray device. Table top tray, the feature will x-ray beams. For each position of the patient, film cassettes of different sizes can be fitted on either side or head and foot end.

Accessories, located on both sides of the table can be easily mounted on stainless steel edge rails Using this system, any accessory, insert or remove the others are not required to remove. The rails are fixed during the entire table is mounted.

The bottom of the table, is designed in the form of " T", so that the doctor, may intervene closely to patient. Table has four antistatic wheel, these are independent wheels rotate in all directions. There is a foot pedal or electro-hydraulic brake system. The table, carrying capacity is 350 kg. When using the table, the weight of the patients, tremors, shaking, interruption, loss of balance, without given any desired position easily. Operating table, in order to ensure the safety of the patient and the doctor protected against static electricity, insulation and grounding protection was necessary.

Optional Features:

- C-arm scope device, for a more efficient use,, one-piece leg section , Head section and leg section can changeable system.
- table can sliding from head to leg 30 cm.
- With manual or electrohydraulic increase /decrease elevator system.
- Electrohydraulic brake system
- Table can rotation 360 degree.
- Accessories for the decline in rail, special wedge assembly
- Visco - elastic cushion set
- 2 Pieces back bearing apparatus
- Safety system for manual using.



STR 2000F model differs from STR 2000 model with a backup battery, extra control panel on base column and sliding function.

Standard Accessories


Polyurethane Cushion set
Arm Support
Shoulder Support
Knee Support
Body Support
Arm Belt
IV Stand
Anesthesia Frame
Waste Tub
Hand Control
Radiographic Benchtop set

Optional Accessories

Single Piece Leg Section
Pin Head, Seating Position belt
Kidney Elevator
Arm Piece
Silicon Head and Body
Support
Wristband
Orthopedic Traction System
"U" Type head
Accessories Table
X-Ray Cassette Holder
Proctology Piece



POSITIONS

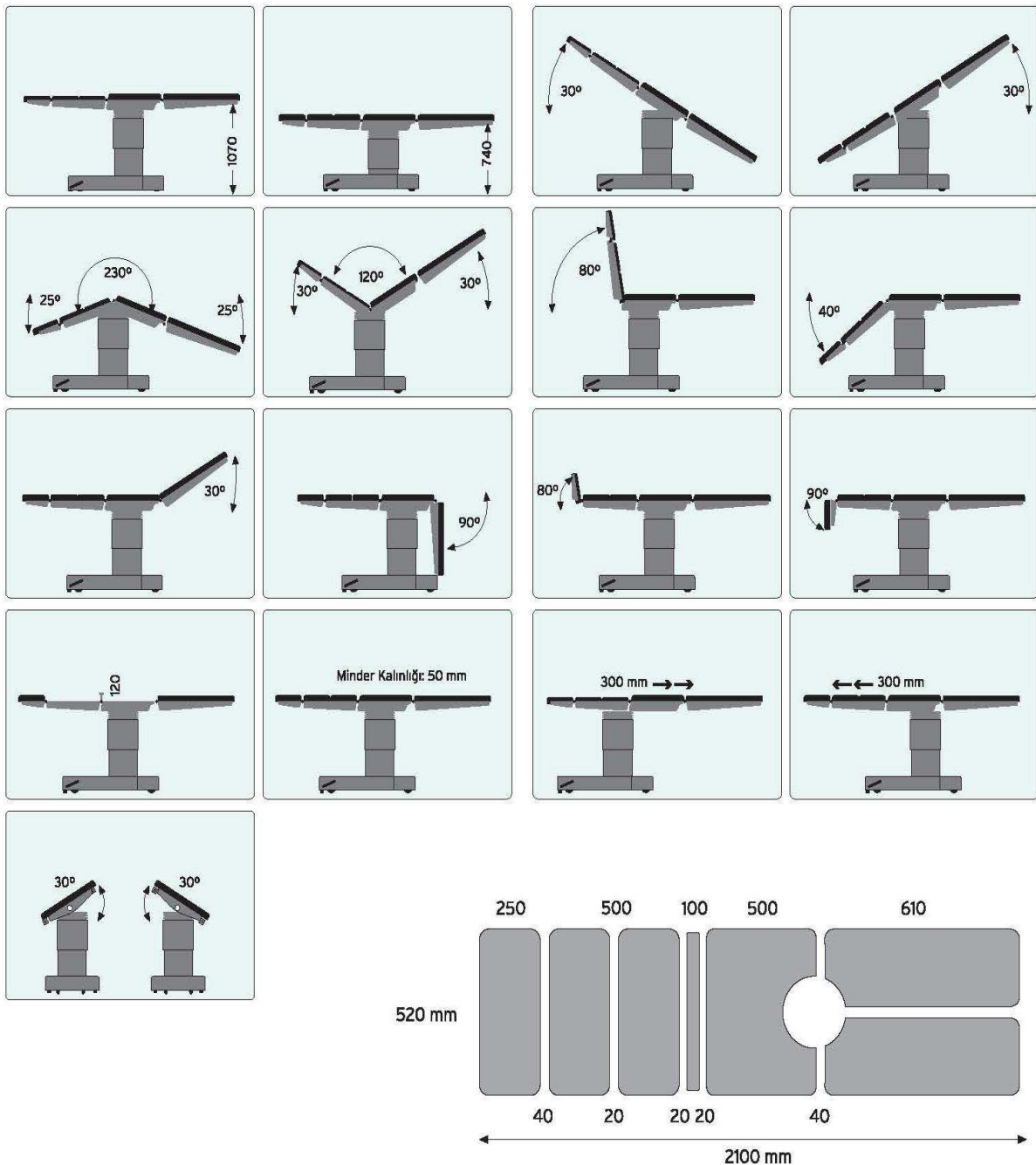

Trendelenburg	30°	Leg section down	90°
Reverse Trendelenburg	30°	Leg section up	30°
Table elevation (min)	740mm	Tilt right – left	25°
Table elevation (max)	1070mm	Fleks (Kidney position)	230°
Back section down	40°	Refleks	120°
Back section up	80°	Kidney Elevator	120 mm
Head section down	90°	Leg section horizontal	180°
Head section up	80°	Forward- back sliding	300 mm



eryiğit

TIBBİ CİHAZLAR AŞ

İvedik Organize Sanayi Bölgesi
Özanadolu Sanayi Sitesi 1453. Sok. No:3
Ostim, Yenimahalle – Ankara
Tel : +90.312.395 57 95
Faks : +90.312.395 57 96
E-posta : info@eryigit.com.tr
Web : www.eryigit.com.tr