

COMEN

Share With The World.



STAR5000C

Maternal/Fetal Monitor

CE 0434





STAR5000C

Maternal/Fetal Monitor

Features


- **Simple operation mode** Simple mode design, 12.1" color TFT touch screen, easy for observation;
- **Superior products performance** Modular design with advanced ARM system, superior module for fetal, ECG, NIBP, convenient for maintenance and upgrade;
- **Excellent system function** Supporting different interfaces like 3-parameter, 9-parameter, big font; FM measurement auto/manual, automatic fetal scoring system;
- **Powerful storage and replaying system** 720 hours all-parameter real time displaying, 30 days trendgraph, 400 groups trend list, 1-hour voice storage, simultaneous review;
- **Alarm system** Intelligent 3-level audio/visual alarm, more accurate analysis;
- **Networking system** USB socket for software upgrade, ethernet socket supporting wire/wireless central monitoring for maximum 128 units;



Folding-up LCD screen, easy for observatio.



112mm wide thermal printer, paper rolls, Z-fold paper are Available.

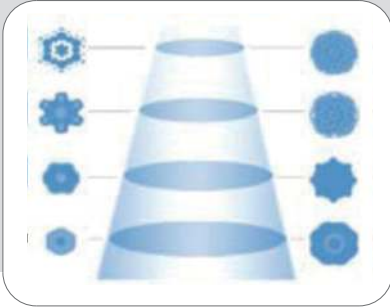

12- crystal pulse wave sensor, for capturing more accurate FHR signal



Hidden handhold, convenient to carry



12.1 inch LCD touch screen, convenient and efficient



External printer is available



Portable wallmount solution



STAR5000C

Maternal/Fetal Monitor

www.szcomen.com

Standard

- FHR
- TOCO
- FM
- ECG/HR
- SPO₂
- NIBP
- Respiration
- Temp
- PR
- Internal Printer

Optional

- Twins
- Neonatal, Pediatric accessories
- Fetal Acoustic Stimulator

Interface



Big Font

Big font display is convenient for monitoring in long distance, long time




Trendgraph

12, 48, 96 hours trendgraph is optional




All-parameter twins monitoring


10 parameters monitoring interface, timely monitoring for both mother and fetus

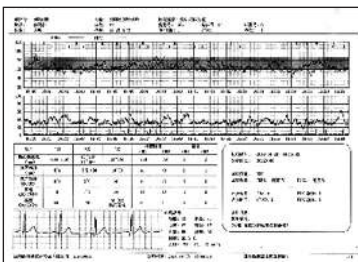

Analysis report print

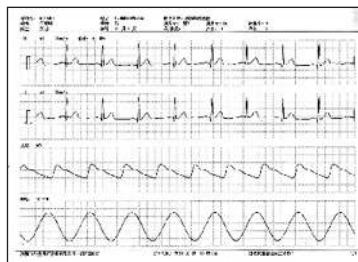

Waveform print

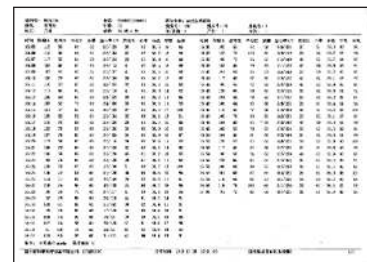

List print

Shenzhen Comen Medical Instruments Co., Ltd.

Add: Floor 7, Block 5, 4th Industrial Park of Nanyou, Nanshan District, Shenzhen City 518052, China

Tel: +86-755-2640 8879 2641 9446 Fax: +86-755-2643 1232

Website: www.szcomen.com E-mail: info@szcomen.com